

Chelating Agent

Column: Luna 5u C18(2) 100A

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4252-E0

Elution Type: Isocratic

Eluent A: 72% Tetrabutylammonium hydroxide, 0.083M, pH 3.9 with phosphoric acid

Pct. A: 72

Eluent B: Acetonitrile


Pct. B: 28

Flow Rate: 1.2 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs. - Variable Wave. @ 254 nm (ambient)

Page 1 of 1

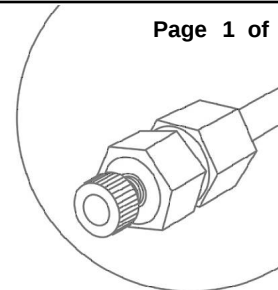

Products used in this Application:

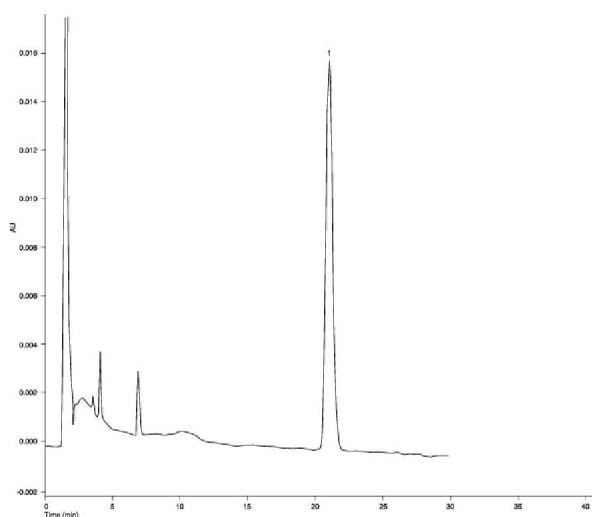LUNA®

ANALYTES:

1 EDTA