



FAMES

Column: Zebron ZB-5MS
Phase: 5% Polysilarylene - 95% Polydimethylsiloxane copolymer
Dimensions: 30 meters x 0.32 mm x 0.25 µm
Order No: 7HM-G010-11
Oven Profile: 100°C for 2 min to 280°C at 8°C/min for 5 min
Carrier Gas: Constant Flow Helium, 1.1 mL/min
Injection: split 20:1 0.2 µL @ 280 °C
Detection: Flame Ionization (FID) (325 °C)

Page 1 of 1

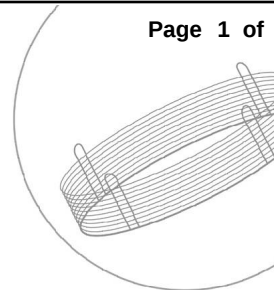

Products used in this Application:

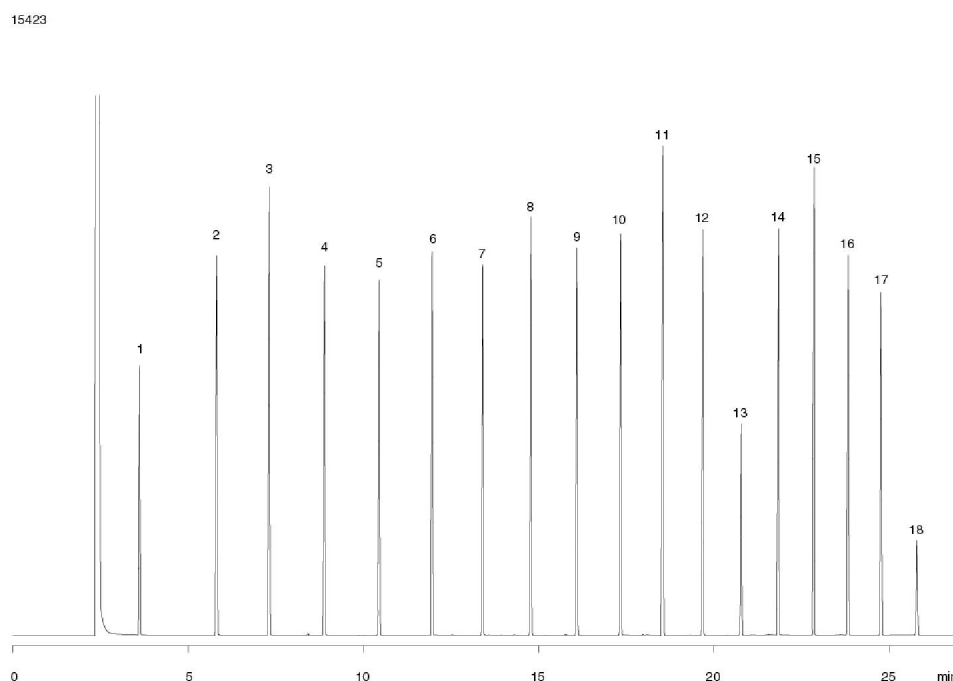Zebron™

ANALYTES:

- | | |
|-------------------------------|--------------------------------|
| 1 Methyl Caproate (C6) | 10 Methyl Palmitate (C16) |
| 2 Methyl Caprylate (C8) | 11 Methyl Heptadecanoate (C17) |
| 3 Methyl Nonoate (C9) | 12 Methyl Stearate (C18) |
| 4 Methyl Caprate (C10) | 13 Methyl Nonadecanoate (C19) |
| 5 Methyl Undecanoate (C11) | 14 Methyl Arachidate (C20) |
| 6 Methyl Laurate (C12) | 15 Methyl Heneicosanoate (C21) |
| 7 Methyl Tridecanoate (C13) | 16 Methyl Behenate (C22) |
| 8 Methyl Myristate (C14) | 17 Methyl Tricosanoate (C23) |
| 9 Methyl Pentadecanoate (C15) | 18 Methyl Lignocerate (C24) |

©Copyright 2000-2007 Phenomenex Inc.

For more information call your Phenomenex Representative at (310) 212-0555

www.phenomenex.com Phenomenex products are available worldwide. Email us at international@phenomenex.com					
Australia tel.: 02-9428-6444 email: info@phenomenex.com.au	Austria 01-319-1301 anfrage@phenomenex.com	Canada (800) 543-3681 info@phenomenex.com	Denmark 4824 8048 dkinfo@phenomenex.com	France 01 30 09 21 10 franceinfo@phenomenex.com	Germany 06021-58830-0 anfrage@phenomenex.com
Ireland tel.: 01 247 5405 email: eireinfo@phenomenex.com	Italy 051 736176 italiainfo@phenomenex.com	New Zealand 09-4780951 info@phenomenex.co.nz	Puerto Rico (800) 541-HPLC info@phenomenex.com	United Kingdom 01625-501367 ukinfo@phenomenex.com	USA (310) 212-0555 info@phenomenex.com