

Diazepam Analysis in Urine

Column: Onyx Monolithic C18

Dimensions: 100 x 4.6 mm ID

Order No: CH0-7643

Elution Type: Gradient

Eluent A: 0.1% formic acid in Water


Eluent B: 0.1% formic acid in Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	4	10	90


Flow Rate: 4 mL/min

Col. Temp.: 30 °C

Detection: Mass Spectrometer (MS) (ambient)

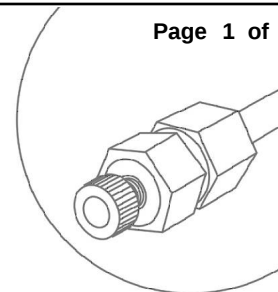

Products used in this Application:

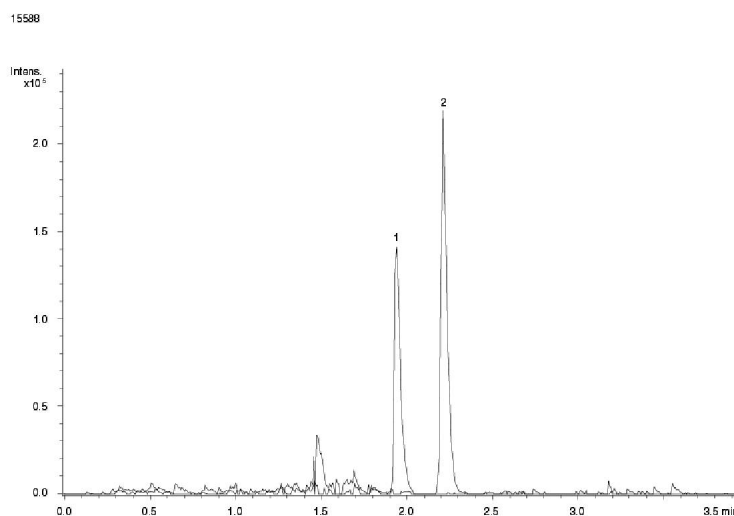

ANALYTES:

- 1 Nordiazepam
- 2 Diazepam