



Goldenseal on Synergi Hydro-RP

Page 1 of 1

Column: Synergi 4u Hydro-RP 80A
Dimensions: 150 x 4.6 mm ID
Order No: 00F-4375-E0
Elution Type: Isocratic
Eluent A: 100mM Dihydrogen Potassium phosphate, pH 4.6 **Pct. A:** 73
Eluent B: Acetonitrile **Pct. B:** 27
Flow Rate: 1.5 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs. - Variable Wave. @ 235 nm (ambient)

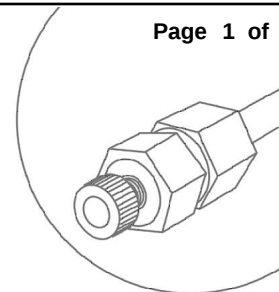

Products used in this Application:

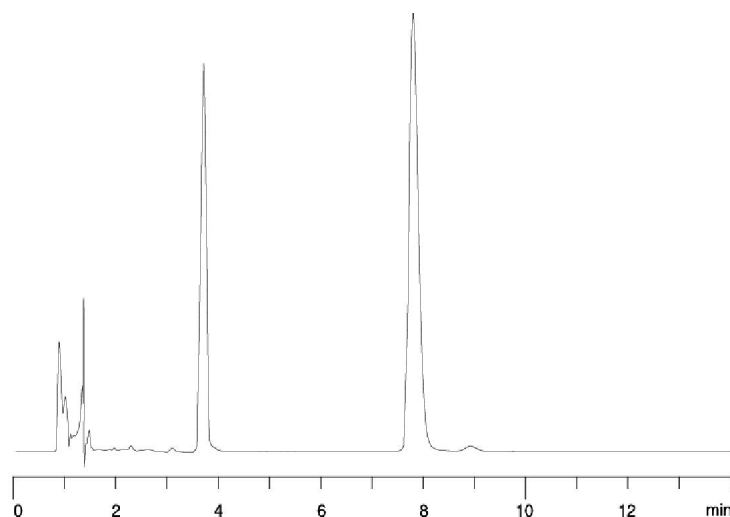synergi[™]

ANALYTES:

- 1 Hydrastatine
- 2 Berberine



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country,
contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.

USA
tel.: (310) 212-0555
fax: (310) 328-7768
email: info@phenomenex.com

Puerto Rico
(800) 541-HPLC
(310) 328-7768
info@phenomenex.com

Canada
(800) 543-3681
(310) 328-7768
info@phenomenex.com

France
01 30 09 21 10
01 30 09 21 11
franceinfo@phenomenex.com

United Kingdom
01625-501367
01625-501796
ukinfo@phenomenex.com

Germany
06021-58830-0
06021-58830-11
anfrage@phenomenex.com

New Zealand
09-4780951
09-4780952
info@phenomenex.co.nz

Australia
02-9428-6444
02-9428-6445
info@phenomenex.com.au