

Ginseng - Luna 5µm C18(2)

Column: Luna 5u C18(2) 100A

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4252-E0

Elution Type: Gradient

Eluent A: Water


Eluent B: Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	85	15
	2	60	20	80


Flow Rate: 1.5 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs. - Variable Wave. @ 205 nm (ambient)

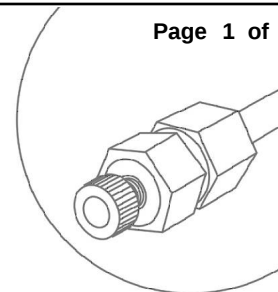

Products used in this Application:

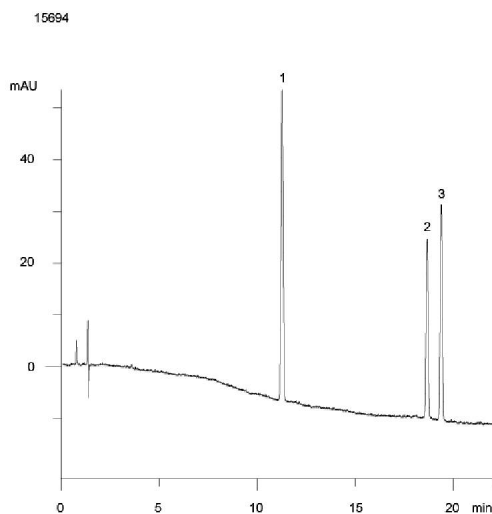LUNA®

ANALYTES:

- 1 Ginsenoside-Re
- 2 Ginsenoside-Rb1
- 3 Ginsenoside-Rc