

Nucleobases

Column: Luna 5u C18 100A

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4041-E0

Elution Type: Isocratic

Eluent A: Acetonitrile

Eluent B: Buffer (10 mM), pH 5.0

Flow Rate: 1 mL/min


Col. Temp.: 25 °C

Detection: UV-Vis Abs. - Variable Wave. @ 212 nm (ambient)

Pct. A: 1


Pct. B: 99

Page 1 of 1

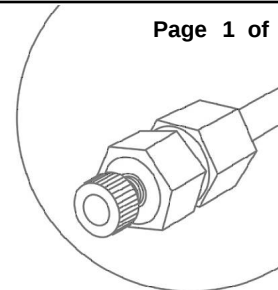

Products used in this Application:

LUNA™

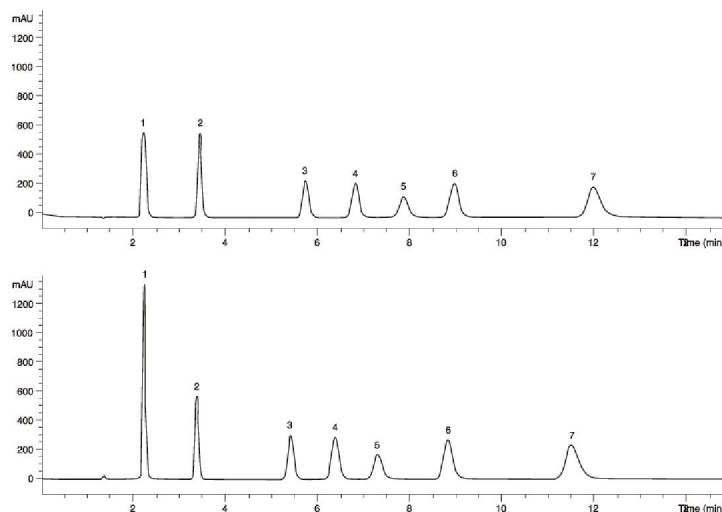

ANALYTES:

- 1 Cytosine
- 2 Uracil
- 3 Adenine
- 4 Guanine
- 5 Xanthine
- 6 Uridine
- 7 Thymine



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.