

Wood saccharides using Rezex RPM

Column: Rezex RPM-Monosaccharide Pb++(8%)

Dimensions: 300 x 7.8 mm ID

Order No: 00H-0135-K0

Elution Type: Isocratic

Eluent A: Water


Pct. A: 100

Flow Rate: 0.6 mL/min

Col. Temp.: 85 °C

Detection: Evaporative Light Scattering (ambient)


Page 1 of 1

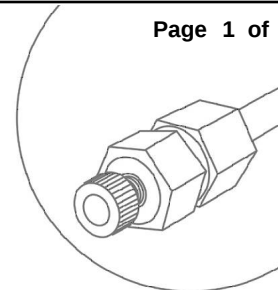

Products used in this Application:

REZEX™

14967

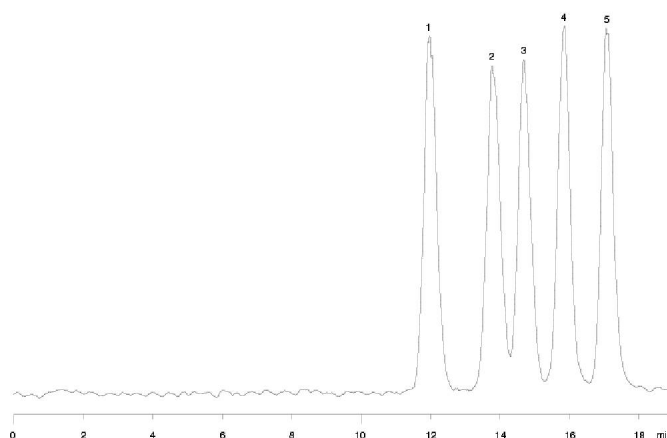

ANALYTES:

- 1 D-Cellulobiose
- 2 α-D-Glucose
- 3 D-(+)-Xylose
- 4 D-(+)-Galactose
- 5 D-(-)-Arabinose