



Blood Alcohols on ZB-WAXplus

Page 1 of 1

Column: Zebron ZB-WAX-Plus
Phase: Polyethylene glycol (PEG)
Dimensions: 30 meters x 0.25 mm x 0.25 µm
Order No: 7HG-G013-11
Oven Profile: 65°C isothermal
Carrier Gas: Constant Flow Helium, 1.5 mL/min
Injection: Split 20:1 50 µL @ 225 °C
Detection: Flame Ionization (FID) (250 °C)

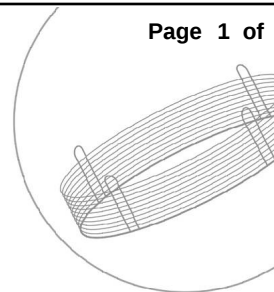

Products used in this Application:

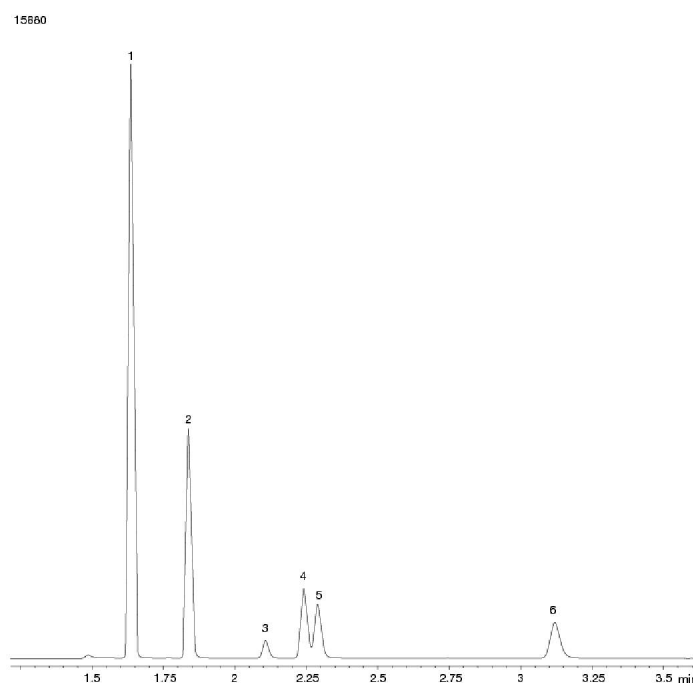Zebron[™]

ANALYTES:

- 1 Acetaldehyde
- 2 Acetone
- 3 Methanol
- 4 Isopropanol
- 5 Ethanol
- 6 n-Propanol

©Copyright 2000-2006 Phenomenex Inc.

For more information call your Phenomenex Representative at (310) 212-0555



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.

phenomenex[®]
...breaking with tradition[™]

USA
tel.: (310) 212-0555
fax: (310) 328-7768
email: info@phenomenex.com

Puerto Rico
(800) 541-HPLC
(310) 328-7768
info@phenomenex.com

Canada
(800) 543-3681
(310) 328-7768
info@phenomenex.com

France
01 30 09 21 10
01 30 09 21 11
franceinfo@phenomenex.com

United Kingdom
01625-501367
01625-501796
ukinfo@phenomenex.com

Germany
06021-58830-0
06021-58830-11
anfrage@phenomenex.com

New Zealand
09-4780951
09-4780952
info@phenomenex.co.nz

Australia
02-9428-6444
02-9428-6445
info@phenomenex.com.au