

Purines and pyrimidines on Luna Amino

Column: Luna 5u NH2 100A

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4378-E0

Elution Type: Isocratic

Eluent A: Acetonitrile

Pct. A: 90

Eluent B: Water


Pct. B: 10

Flow Rate: 2 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs. - Variable Wave. @ 254 nm (ambient)

Page 1 of 1

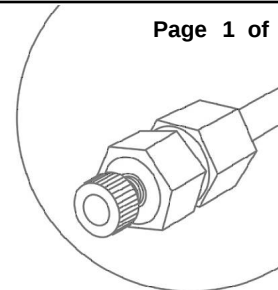

Products used in this Application:

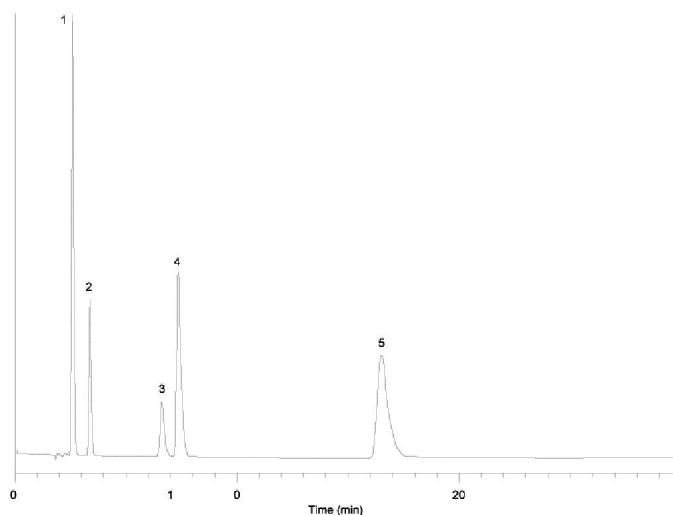LUNA[™]

ANALYTES:

- 1 Thymine
- 2 Uracil
- 3 Cytosine
- 4 Adenine
- 5 Guanine



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.