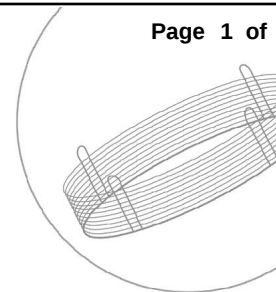



**Abused Inhalants ZB-624 30x0.32x1.8**

Page 1 of 1

**Column:** Zebron ZB-624  
**Phase:** 6%-cyanopropylphenyl-94%-methyl polysiloxane  
**Dimensions:** 30 meters x 0.32 mm x 1.8 µm  
**Order No:** 7HM-G005-31  
**Oven Profile:** 40°C for 4 min to 120°C at 10°C/min for 4 min.  
**Carrier Gas:** Constant Flow Helium, 3.5 mL/min  
**Injection:** Split 20:1 100 µL @ 125 °C  
**Detection:** Flame Ionization (FID) (250 °C)

**Products used in this Application:****Zebron™****Analyst Note:** Headspace analysis

14860

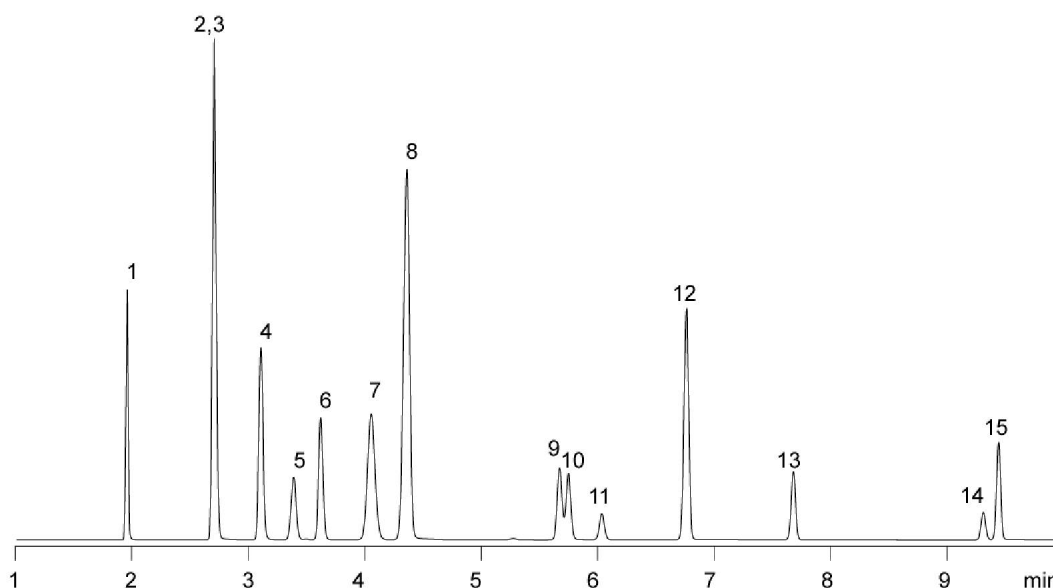**ANALYTES:**

- |                                  |                                  |
|----------------------------------|----------------------------------|
| 1 Methanol                       | 10 Ethyl acetate                 |
| 2 Ethanol                        | 11 Chloroform                    |
| 3 Ether                          | 12 Benzene                       |
| 4 Acetone                        | 13 Trichloroethylene             |
| 5 Isopropanol                    | 14 Methyl Isobutyl Ketone (MIBK) |
| 6 Methylene chloride             | 15 Toluene                       |
| 7 Methyl tert butyl ether (MTBE) |                                  |
| 8 Hexane                         |                                  |
| 9 Methyl ethyl ketone (MEK)      |                                  |

©Copyright 2000-2008 Phenomenex Inc.

For more information call your Phenomenex Representative at (310) 212-0555

|                                                                                                                                                                   |                                                                                                        |                                                                                                      |                                                                                                      |                                                                                                             |                                                                                                       |
|-------------------------------------------------------------------------------------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------------------------|------------------------------------------------------------------------------------------------------|------------------------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------------------------|
| <b>www.phenomenex.com</b> Phenomenex products are available worldwide. Email us at <a href="mailto:international@phenomenex.com">international@phenomenex.com</a> |                                                                                                        |                                                                                                      |                                                                                                      |                                                                                                             |                                                                                                       |
| <b>Australia</b><br>tel.: 02-9428-6444<br>email: <a href="mailto:info@phenomenex.com.au">info@phenomenex.com.au</a>                                               | <b>Austria</b><br>01-319-1301<br><a href="mailto:anfrage@phenomenex.com">anfrage@phenomenex.com</a>    | <b>Canada</b><br>(800) 543-3681<br><a href="mailto:info@phenomenex.com">info@phenomenex.com</a>      | <b>Denmark</b><br>4824 8048<br><a href="mailto:dkinfo@phenomenex.com">dkinfo@phenomenex.com</a>      | <b>France</b><br>01 30 09 21 10<br><a href="mailto:franceinfo@phenomenex.com">franceinfo@phenomenex.com</a> | <b>Germany</b><br>06021-58830-0<br><a href="mailto:anfrage@phenomenex.com">anfrage@phenomenex.com</a> |
| <b>Ireland</b><br>tel.: 01 247 5405<br>email: <a href="mailto:eireinfo@phenomenex.com">eireinfo@phenomenex.com</a>                                                | <b>Italy</b><br>051 736176<br><a href="mailto:italiainfo@phenomenex.com">italiainfo@phenomenex.com</a> | <b>New Zealand</b><br>09-4780951<br><a href="mailto:info@phenomenex.co.nz">info@phenomenex.co.nz</a> | <b>Puerto Rico</b><br>(800) 541-HPLC<br><a href="mailto:info@phenomenex.com">info@phenomenex.com</a> | <b>United Kingdom</b><br>01625-501367<br><a href="mailto:ukinfo@phenomenex.com">ukinfo@phenomenex.com</a>   | <b>USA</b><br>(310) 212-0555<br><a href="mailto:info@phenomenex.com">info@phenomenex.com</a>          |