



ASTM Method D5135-95: Solvent Analysis

Page 1 of 1

Column: Zebron ZB-WAX
Phase: Polyethylene glycol (PEG)
Dimensions: 30 meters x 0.25 mm x 0.25 µm
Order No: 7HG-G007-11
Oven Profile: 35°C hold for 3min. to 95°C at 20°C/min. hold for 6min.
Carrier Gas: Constant Pressure Helium, 1.6 mL/min
Injection: Split 100:1 0.5 µL @ 250 °C
Detection: Flame Ionization (FID) (260 °C)

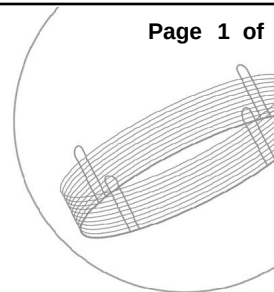

Products used in this Application:

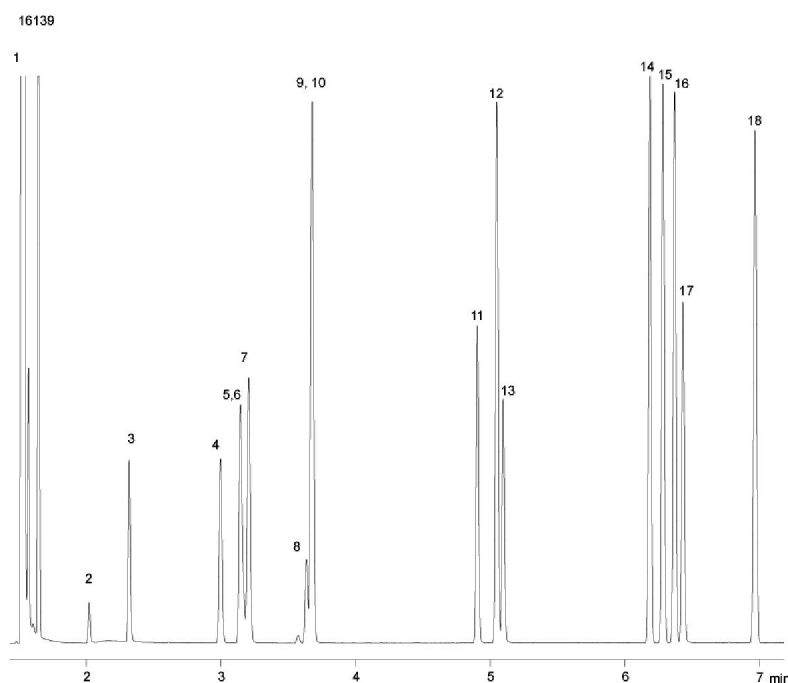Zebron™

ANALYTES:

- | | |
|-----------------------|-----------------|
| 1 Pentane | 10 Benzene |
| 2 Methyl Formate | 11 2-Butanol |
| 3 Acetone | 12 Toluene |
| 4 Ethyl Acetate | 13 Propanol |
| 5 Methyl ethyl ketone | 14 Ethylbenzene |
| 6 Methanol | 15 p-Xylene |
| 7 2-Methyl-2-propanol | 16 m-Xylene |
| 8 Methylene chloride | 17 1-Butanol |
| 9 Ethanol | 18 o-Xylene |

©Copyright 2000-2008 Phenomenex Inc.

For more information call your Phenomenex Representative at (310) 212-0555

www.phenomenex.com Phenomenex products are available worldwide. Email us at international@phenomenex.com					
Australia tel.: 02-9428-6444 email: info@phenomenex.com.au	Austria 01-319-1301 anfrage@phenomenex.com	Canada (800) 543-3681 info@phenomenex.com	Denmark 4824 8048 dkinfo@phenomenex.com	France 01 30 09 21 10 franceinfo@phenomenex.com	Germany 06021-58830-0 anfrage@phenomenex.com
Ireland tel.: 01 247 5405 email: eireinfo@phenomenex.com	Italy 051 736176 italiainfo@phenomenex.com	New Zealand 09-4780951 info@phenomenex.co.nz	Puerto Rico (800) 541-HPLC info@phenomenex.com	United Kingdom 01625-501367 ukinfo@phenomenex.com	USA (310) 212-0555 info@phenomenex.com