

Polyethylene Glycol 330

Column: Phenogel 5u 50A

Dimensions: 300 x 7.8 mm ID

Order No: 00H-0441-K0

Elution Type: Isocratic


Eluent A: Tetrahydrofuran

Pct. A: 100

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: Differential refractometer (ambient)

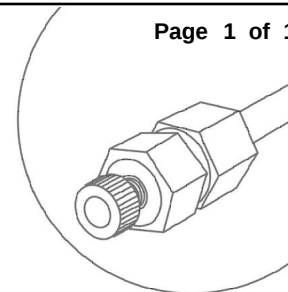


Page 1 of 1

Products used in this Application:

Phenogel™

Analyst Note: 3 Phenogel 5μ columns (300 x 7.8mm ID) placed in series as follows: 50Å, 100Å, 500Å.

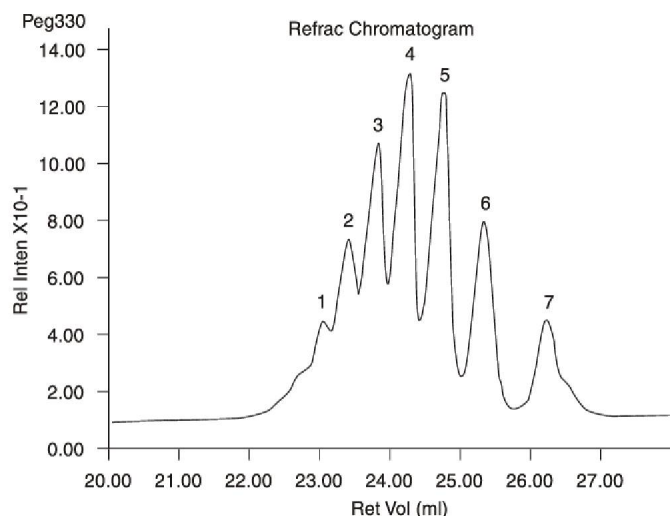

ANALYTES:

- 1 dp 7 546 MW
- 2 dp 6 458 MW
- 3 dp 5 370 MW
- 4 dp 4 282 MW
- 5 dp 3 194 MW
- 6 dp 2 106 MW
- 7 dp 1 62 MW