


## Oligosaccharides using Rezex RSO


**Column:** Rezex RSO-Oligosaccharide Ag+ 4%  
**Dimensions:** 200 x 10 mm ID  
**Order No:** 00P-0133-N0  
**Elution Type:** Isocratic  
**Eluent A:** Water  
**Flow Rate:** 0.3 mL/min  
**Col. Temp.:** 75 °C  
**Detection:** Evaporative Light Scattering (ambient)

Pct. A: 100

Page 1 of 1

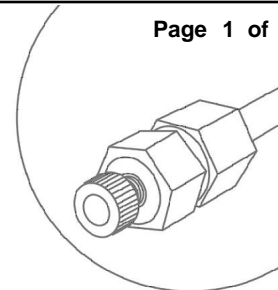

Products used in this Application:

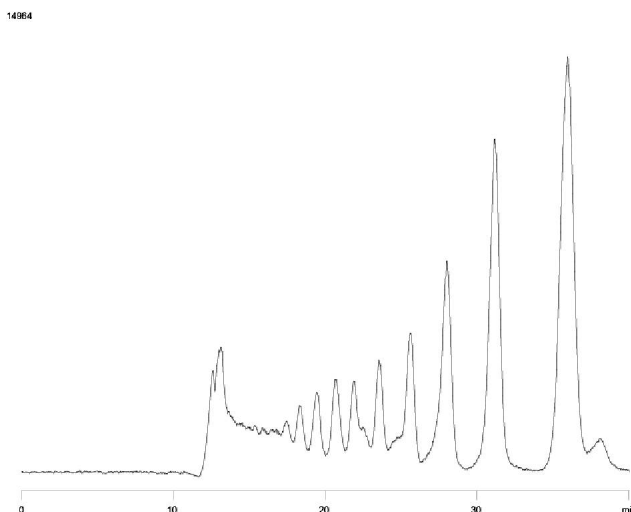**REZEX™**

### ANALYTES:

- 1 Malto-Oligosaccharides Dp1-Dp14