

Saccharide analysis using Rezex RSO

Column: Rezex RSO-Oligosaccharide Ag+ 4%

Dimensions: 300 x 7.8 mm ID

Order No: 00H-0133-K0

Elution Type: Isocratic

Eluent A: Water


Flow Rate: 0.6 mL/min

Col. Temp.: 75 °C

Detection: Evaporative Light Scattering (ambient)


Pct. A: 100

Page 1 of 1

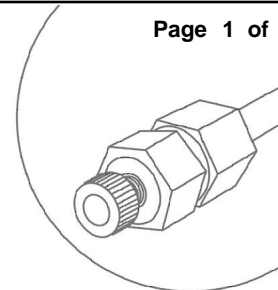

Products used in this Application:

REZEX[™]

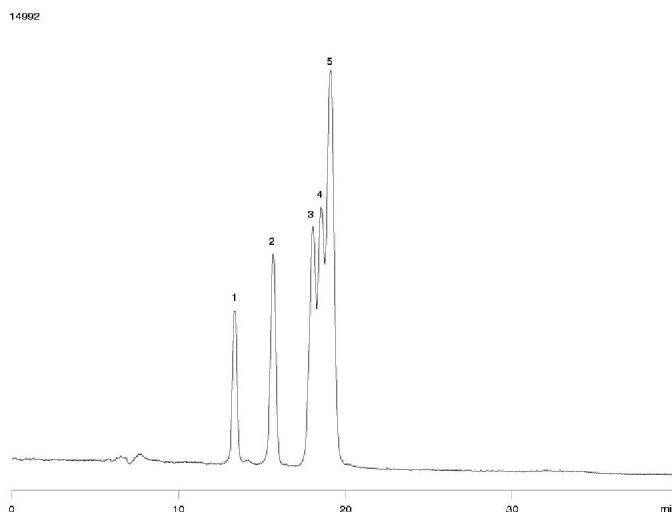

ANALYTES:

- 1 Melezitose
- 2 Maltose
- 3 Glucose
- 4 Mannose
- 5 Fructose
- 6 Ribitol