


Acidic compounds on Synergi Polar-RP

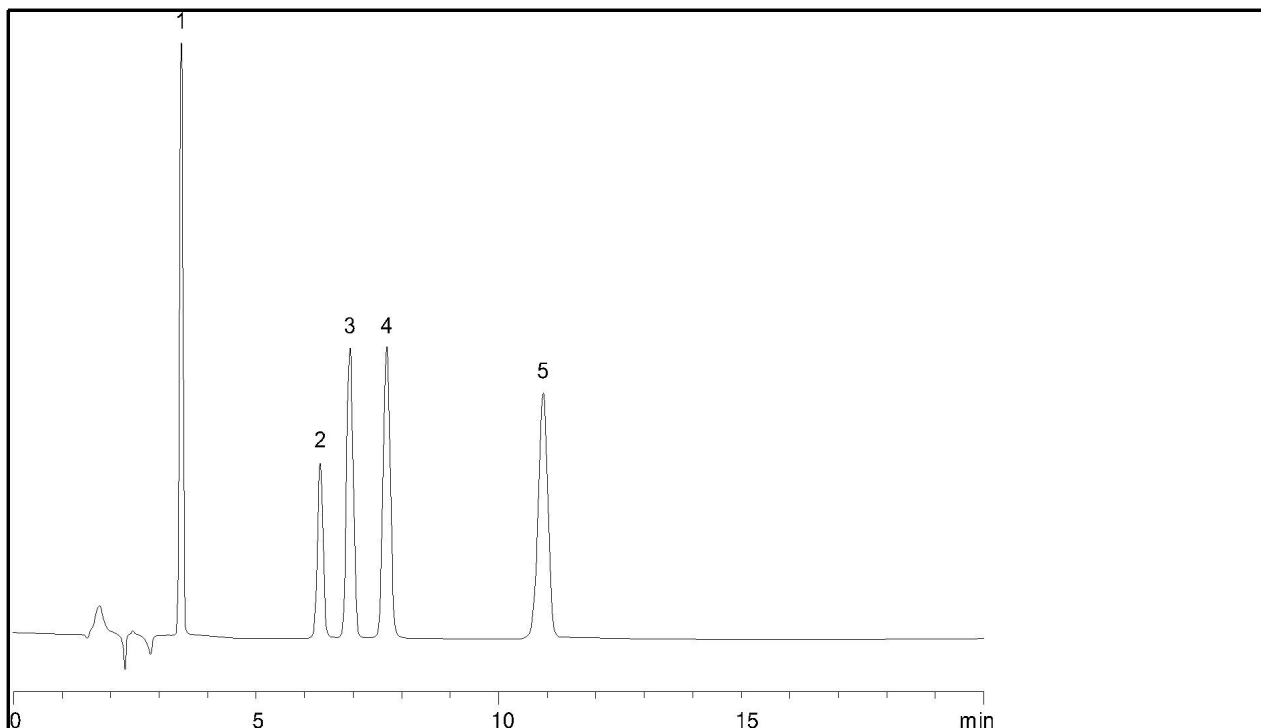

Column: Synergi 4 μ m C6H5 / Phenyl Polar-RP

Dimensions: 150 x 4.6 mm ID

Order No:

Elution Type: Isocratic

Eluent A: 20mM Potassium phosphate, pH 2.5

Pct. A: 75

Eluent B: Acetonitrile

Pct. B: 25

Flow Rate: 1.5 mL/min

Temperature: 30 °C

Detection: UV-Vis Abs. - Variable Wave. @ 254 nm

SAMPLE ANALYTE

- 1 *p*-Hydroxybenzoic acid
- 2 Sorbic acid
- 3 Benzoic acid
- 4 Salicylic acid
- 5 *p*-Toluic acid

SPE:

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 14273