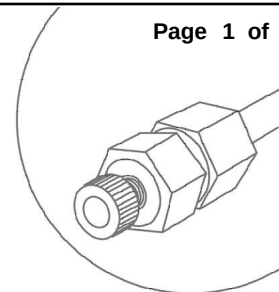
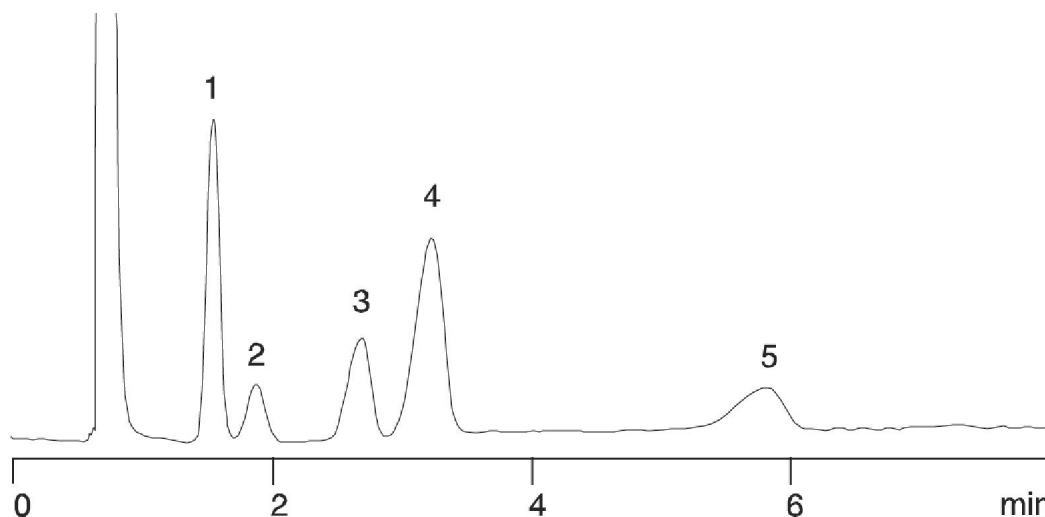


Simple Sugars (#6) on Luna Amino

Column: Luna 5u NH2 100A**Dimensions:** 250 x 4.6 mm ID**Order No:** 00G-4378-E0**Elution Type:** Isocratic**Eluent A:** Acetonitrile**Pct. A:** 80**Eluent B:** Water**Pct. B:** 20**Flow Rate:** 3 mL/min**Col. Temp.:** 40 °C**Detection:** Refractive Index (RI) (ambient)

Page 1 of 1

**Products used in this Application:****LUNA®****ANALYTES:**

- 1 Fructose
- 2 Dextrose
- 3 Sucrose
- 4 Maltose
- 5 Maltotriose

www.phenomenex.comPhenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.