


Sugar Alcohols on Luna Amino

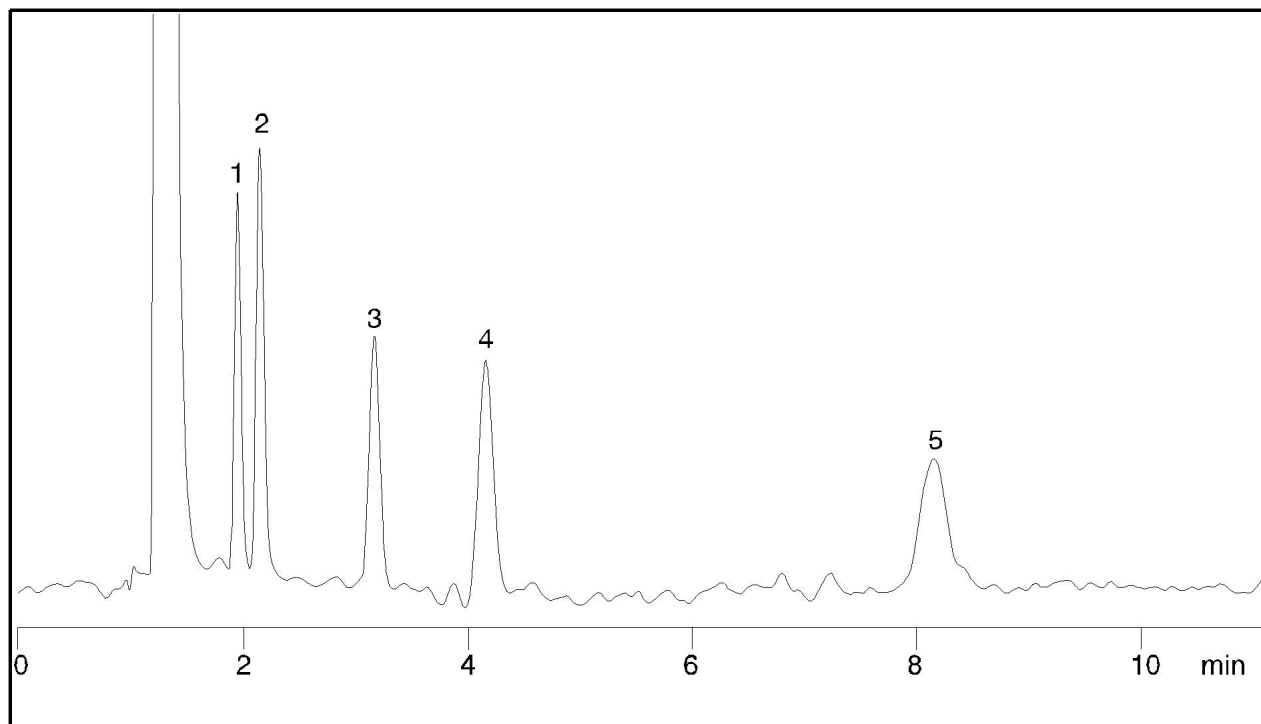

Column: Luna 5 μ m NH₂

Dimensions: 250 x 4.6 mm ID

Order No:

Elution Type: Isocratic

Eluent A: Acetonitrile

Pct. A: 85

Eluent B: Water

Pct. B: 15

Flow Rate: 3 mL/min

Temperature: 40 °C

Detection: Refractive Index @

SAMPLE ANALYTE

- 1 Glycerol
- 2 Pentaerythritol
- 3 Xylitol
- 4 Sorbitol
- 5 Inositol

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. 14289
No.