


Steroids

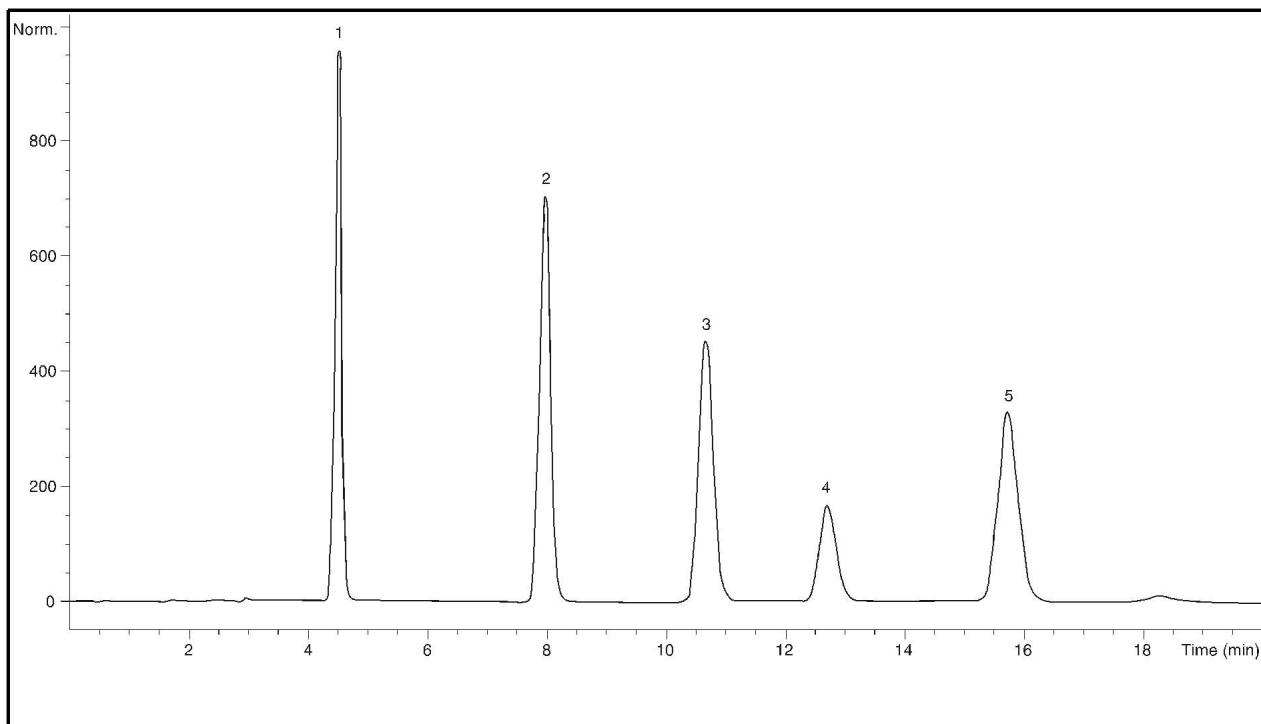

Column: Luna 5 μ m C18 / ODS Original

Dimensions: 150 x 4.6 mm ID

Order No:

Elution Type: Isocratic

Eluent A: 0.1% H₃PO₄

Pct. A: 54

Eluent B: Acetonitrile

Pct. B: 35

Eluent C: Methanol

Pct. C: 11

Flow Rate: 1 mL/min

Temperature: 27 °C

Detection: UV-Vis Abs. - Variable Wave. @ 254 nm

SAMPLE ANALYTE

- 1 Hydrocortisone
- 2 Corticosterone
- 3 11-Alpha-Hydroxyprogesterone
- 4 Cortisone acetate
- 5 11-Ketoprogesterone

SPE:

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 1185/ALG/34