

## Naproxen Tablets - USP Method

Column: Synergi 4u Fusion-RP 80A

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4424-E0

Elution Type: Isocratic

Eluent A: Water

Eluent B: Methanol

Eluent C: Glacial Acetic Acid

Flow Rate: 1.2 mL/min

Col. Temp.: 25 °C


Detection: UV-Vis Abs. - Variable Wave. @ 254 nm (25 °C)

Pct. A: 49


Pct. B: 50

Pct. C: 1

Page 1 of 1

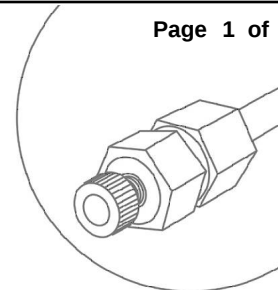

Products used in this Application:

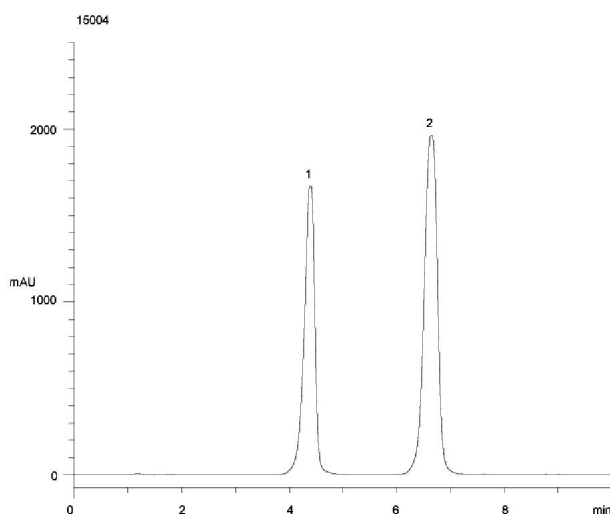synergi<sup>™</sup>

## ANALYTES:

- 1 Naproxen
- 2 Butyrophenone



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.

**USA**  
tel.: (310) 212-0555  
fax: (310) 328-7768  
email: info@phenomenex.com

**Puerto Rico**  
(800) 541-HPLC  
(310) 328-7768  
info@phenomenex.com

**Canada**  
(800) 543-3681  
(310) 328-7768  
info@phenomenex.com

**France**  
01 30 09 21 10  
01 30 09 21 11  
franceinfo@phenomenex.com

**United Kingdom**  
01625-501367  
01625-501796  
ukinfo@phenomenex.com

**Germany**  
06021-58830-0  
06021-58830-11  
anfrage@phenomenex.com

**New Zealand**  
09-4780951  
09-4780952  
info@phenomenex.co.nz

**Australia**  
02-9428-6444  
02-9428-6445  
info@phenomenex.com.au