


Urea-Glycerol on Luna Amino

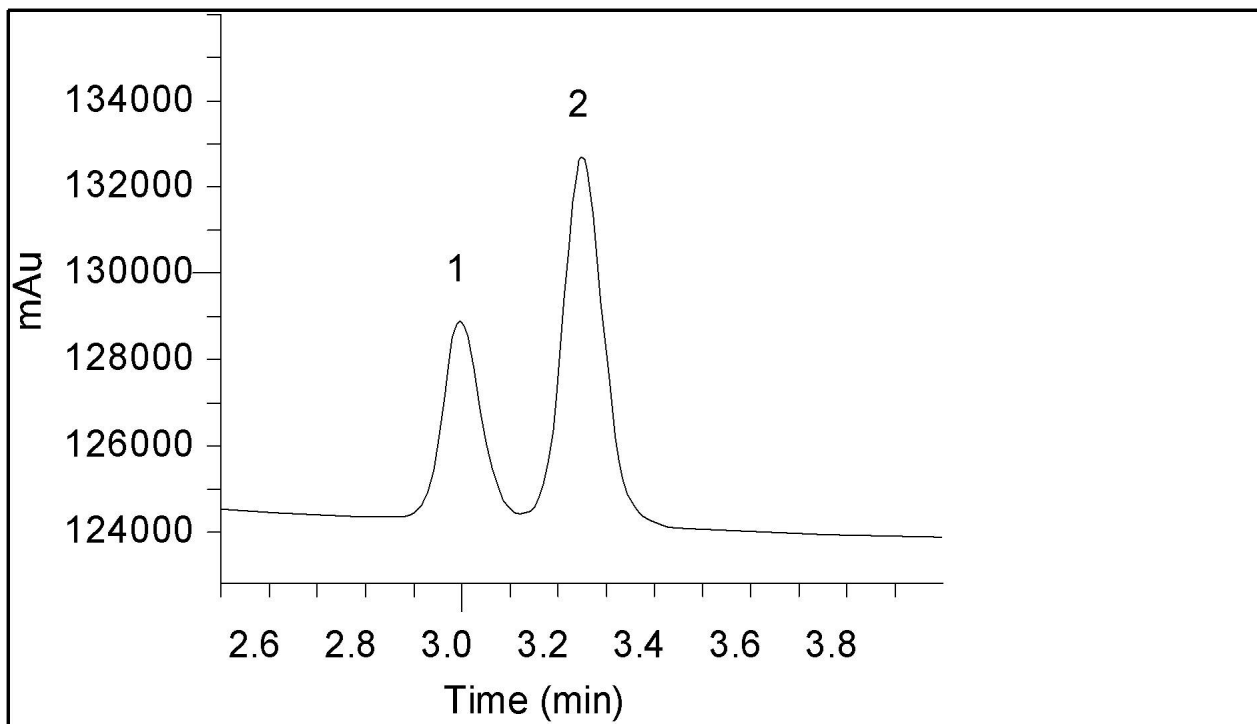

Column: Luna 5 μ m NH₂

Dimensions: 250 x 4.6 mm ID

Order No:

Elution Type: Isocratic

Eluent A: 90:10 Acetonitrile:Water

Pct. A: 100

Flow Rate: 2 mL/min

Temperature: 22 °C

Detection: Refractive Index @

SAMPLE ANALYTE

1 Urea

2 Glycerol

SPE:

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 14362