


Tricyclic Antidepressants at pH 11 in Gradient RP HPLC



15017

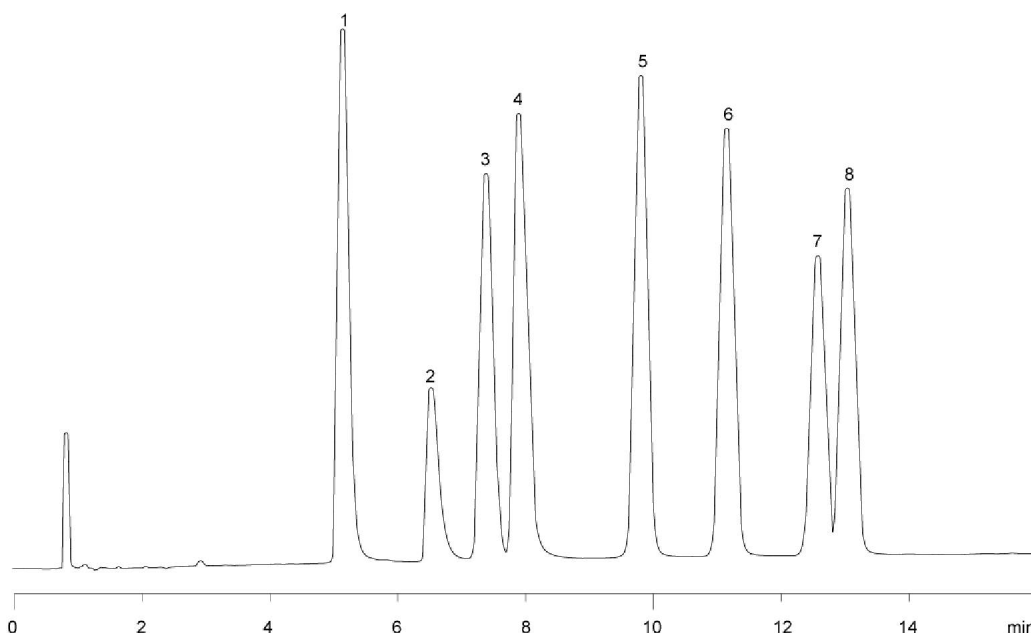

Column: Gemini 5 μ m C18 / ODS C18
Dimensions: 150 x 4.6 mm ID
Order No: 00F-4435-E0
Elution Type: Gradient
Eluent A: 20mM TEA, pH 11.0
Eluent B: ACN
Eluent C: MeOH

Gradient Profile:

Step No.	Step Time	Pct A	Pct B	Pct C
1	0	35	25	40
2	15	20	25	55

Flow Rate: 1.5 mL/min
Temperature: 22 °C
Detection: UV-Vis Abs. - Variable Wave. @ 254 nm
SPE:

SAMPLE ANALYTES:

- 1 Nordoxepine
- 2 Protriptyline
- 3 Doxepine
- 4 Nortriptyline
- 5 Imipramine
- 6 Amitriptyline
- 7 Clomipramine
- 8 Trimipramine

For more information call your Phenomenex Representative at

©Copyright 2000-2005 Phenomenex Inc.

Ref. No.: 15017

Page 1 of 1